




NEW SELF-TEST GENERATION

- Detects omicron variant and shows no cross-reactivity
- High quality packaging protects against external influences
- Sealed buffer liquid and test cassette guarantee long shelf life
- BfArM-listed (AT1276/21) and PEI-evaluated (German Health Authorities)
- With CE certification (CE 1011) - Approved for self-testing in EU Complies with EN ISO 13485:2016
- High performance confirmed in most recent evaluation study by RKI (Robert Koch Institute) and PEI in Dec 2021*
**In direct comparison to other positively evaluated rapid tests*
- Easy handling due to prefilled buffer liquid and sample collection through short nose (approx. 1.5 cm):

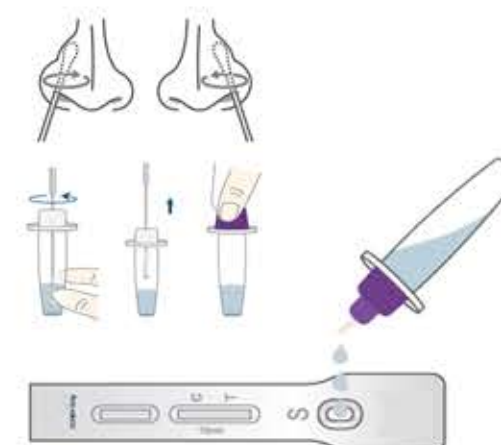

References

- 2,500+ pharmacies (DE and AT)
- Critical infrastructure (DE) / International airports (DE) / DAX corporations
- Leading drugstore chain with 2000+ stores (DE)
- Professional version in long-term use in a wide range of medical facilities, laboratories and testing sites in 10+ European countries

JOINSTAR CE SELF-TEST (COLLOIDAL GOLD)

INSTRUCTIONS FOR USE IN DE | EN | FR | IT | ES | PT

The COVID-19 Antigen Assay (colloidal gold) is a qualitative immunoassay-based test for the detection of SARS-CoV-2 nucleocapsid protein antigens in swab specimens. During this assay, the sample reacts with an anti-SARS-CoV-2 antibody conjugate on the pad coated with particles and then the mixture migrates up the membrane chromatogram by capillary action and reacts with a SARSCoV-2 antibody in the detection zone. If the sample contains SARSCoV-2, a colored line appears in the test line area indicating a positive result. If the sample does not contain SARSCoV-2, no colored line appears in this area, indicating a negative result. As a process control, a colored line always appears in the control line area.

Technical Data

Joinstar COVID-19 Antigen Rapid Test (Colloidal Gold) - CE-Certified Self-Test (Single)

- Sensitivity: 96,1% / Specificity: 99,2%
- Shelf life: 2 years / Storage: Between +2°C to 30°C

Contents and Logistics

All necessary components included.

Scope of delivery for single test (= 1 softpack):

- Coronavirus antigen test cassette
- Sample tube with cap and extraction buffer
- Sterile swab
- Disposable bag for disposal
- Instructions for use in English, French, German, Italian, Spanish and Portuguese

Packed in softpacks: 1 test/kit; 400 kits per carton
Carton content: 400 tests (H: 29 cm, W: 39 cm, L: 59 cm)
Weight single item: 0,026 KG / Carton: 10,4 KG
Euro pallet: 9600 tests
4 cartons/layer; 6 layers high; height: ca 185 cm
32 pallets per truck = 307.200 tests per truck

Our Guarantee

1. Complete documentation with LOT numbers for all tests
2. Temperature-controlled storage (Tests permanently lose functionality when exposed to temperatures below 2 degrees Celsius)
3. Strict quality control for each batch / Original products directly from the manufacturer
4. 30 years of diagnostic know-how with professional after-sales service

CHATEAU NATAL GmbH

Linke Wienzeile 118/19
AT-1060 Vienna

FN 54940 w
UID: ATU38968104

+43 (0) 720 145 296 (INT)
+49 (0) 305 201 4206 (DE)

info@chateau-natal.com
www.chateau-natal.com